



Targeted Refrigeration

The new generation of ultra-low temperature refrigeration system can ensure fast refrigeration and was awarded with the Second Prize under State Technological Invention Award.

Five Magic Weapons for Energy Saving

Imported high-efficiency compressor + separated evaporator + composite heat exchanger + optimized fractional condensation & separation system + high-efficiency thermal insulation system can save energy and reduce power consumption of the freezer by 40%.

Low Noise Design

Low-noise direct cooling circuit + suspension frame & sound absorbing compartment & low-noise fan can reduce noise generated by the entire freezer to the largest extent.



The six sides of the cabinet are made from high-efficiency VIP vacuum thermal insulation plate, and the thermal insulation design of the inner door made from foaming material and the outer door system was awarded with multiple patents, which can improve the thermal insulation performance of the freezer.

Optional Accessories(WiFi/Bluetooth printer/Magnetic lock/Swipe card/Fingerprint)

Connect WiFi to upload data to the cloud space, monitor the device, check the status of the freezer faster, and ensure the safety of the samples;

Optional Bluetooth printer, real-time export temperature data, temperature graph, setting operation and other data, to understand the freezer operating status;

Optional functions such as card unlocking, fingerprint unlocking, and facial recognition unlocking are available to facilitate quick access to samples.







7" HD Intelligent Screen Control System

- ·7-inch HD intelligent display touch screen, userfriendly interface and more accurate temperature control, with Bluetooth, WiFi and other functions, more convenient operation and more comprehensive functions;
- · Can add multiple user accounts to log in and record usage/operation records;
- The high-precision microcomputer control system and platinum resistor sensors enable users to set temperature inside the cabinet within a range from -40 to -86°C;
- Optional functions such as card unlocking, fingerprint unlocking, and facial recognition unlocking are available to facilitate quick access to samples.



Security System

- · The perfect audible & visual alarm system: high and low temperature, sensor failure, power failure, low battery, door ajar, main board communication error etc;
- \cdot Equipped with alarm self-recovery function, it reminds users of unhandled alarm problems in freezer many times;
- The compressor start delay and stopping interval protection can ensure reliable operation;
- · Both the touch screen controller and key-board controller have password protection which can prevent any adjustment of operation without permission.



Refrigeration System

- ·The imported high-efficiency compressor and EBM fan are energy-saving and highly efficient;
- The large finned condenser with a space between fins≤2mm, providing effective heat dissipation.
- · DW-HL398HC with twin-compressor, if one is damaged, the other one could keep the temperature at -70°C stably;





Supplies for Cryopreservation

Freezer racks/boxes for cryopreservation are optional.

Remote alarm contact, RS485

Equipped with RS485 and remote alarm interface to realize remote observation of freezer operation status, real-time reception of product operation information, and ensure the storage safety of samples.



Data Storage

The temperature data recorded by the system can be exported, the operating data of the system settings, the alarm data, and the temperature data can be recorded in 2 minutes, which can effectively record the running status of the product for you for a long time.



Thermal Insulation System

- The two-layer heat insulating foamed door with inner and outer door seal and the insulation design of the outer door system with multiple patents can prevent loss of refrigerating capacity in an effective way;
- The six sides of the cabinet are made from high-performance VIP vacuum insulation material, improving thermal insulation performance to a large extent.




Human-oriented

- ·The new designed assisting handle and vacuum release port let the freezer be operated single-handed to open and close the outer door;
- · The newly added file box makes recording easier and more convenient.
- ·The liner made from Galvanized steel sheet/ Stainless steel for medical use is low-temperature tolerant and corrosion-resistant, which has a long service life and is easy to clean;
- · The universal casters and leveling feet design are more convenient for movement and fixation.


Scope of Application

Suitable for use in blood banks, hospitals, health and disease prevention systems, research institutions, colleges & universities, the electronic industry, biological engineering, laboratories in colleges & universities, military enterprises, deep-sea fishing companies, etc.

Performance Data / Cooling Curve

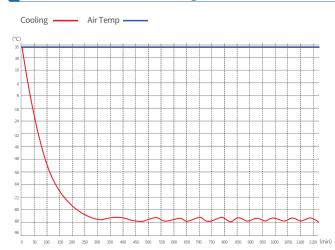

External Dimension

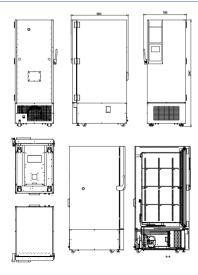

adal ———————————————————————————————————	DW-HL398HC
odel	Upright
abinet Type	398
apacity(L)	440*696*1266
ternal Size(W*D*H)mm	
cternal Size(W*D*H)mm	785*1041*1947
ackage Size(W*D*H)mm	890*1165*2158 247/324
N/GW(Kgs)	241/324
erformance	40, 0000
emperature Range	-40~-86°C 16~32°C
mbient Temperature	
ooling Performance	-86°C
imate Class ontroller	N
	Microprocessor
isplay	HD Intelligent Touch Screen
frigeration	
ompressor	2pcs
poling method	Direct Cooling
efrost Mode	Manual
efrigerant	HC
sulation Thickness(mm)mm	130
onstruction	
kternal Material	Cold rolled steel sheet
ner Material	Galvanized steel sheet/Stainless steel
ielves	3(stainless steel)
oor Lock with key	Yes
adlock	Yes
ccess test port	2pcs. Ø 25 mm
asters	4+(2 leveling feet)
ata Logging/Interval/Recording Time	USB/Record every 5 minutes / 10 years
ackup Battery	Yes
FI	Yes
rm	
mperature	High/Low temperature, High ambient temperature
ectrical	Power failure,Low battery
rstem	Sensor failure, Main board communication error, Condenser overheating alarm, Door ajar
ectrical	
ower Supply(V/HZ)	220-240~/50
ated Current(A)	8.69
ccessories	
tandard	Remote alarm contact,RS485
ptions	Chart recorder, CO ₂ backup system
	L DALL THEOTOPE L LIS DACKIID SVSTAM

JoJo Life Science UG (haftungsbeschränkt) - Biberstraße 32 - 89537 Giengen Tel. 07322-9111329 - Mail: info@jojo-ls.de - Web: www.jojo-ls.de

^{*}The model, parameters and performance specified in this brochure may be changed without prior notice because of product upgrading.

*There may be differences between the product images shown in this brochure and the actual products. When you are buying any product, please check the actual product.